

Newborn Screening

I.	Introduction.....	IV - 1
II.	Test Methods.....	IV - 2
III.	Specimen Collection	IV - 3
	A. Collection Form	IV - 3
	B. Important Information	IV - 3
	C. Instructions for filling out the form.....	IV - 5
	1. Infant's Information.....	IV - 6
	2. Physician or Provider Information.....	IV - 9
	3. Mother's Information.....	IV - 9
	4. Hearing Screen	IV - 10
	5. Date Rec'd/Lab No. (White area at bottom of form)	IV - 11
	D. Filter Paper Collection Procedure	IV - 11
IV.	Shipment of Specimens.....	IV - 13
V.	Reporting and Interpretation.....	IV - 13
VI.	Reference Ranges	IV - 14
VII.	Specimen Criteria	IV - 15
	Chart IV-1, Common Causes of Unsatisfactory Newborn Screening Specimens.....	IV - 16
VIII.	Confirmation Procedures	IV - 17
IX.	Material for Parents	IV - 18
X.	Infant Discharged before Specimen Was Obtained.....	IV - 18
XI.	Refusal to Have Newborn Testing Performed	IV - 18
XII.	Death Notice.....	IV - 18
XIII.	Specimen Collection for the Hemoglobin Electrophoresis Meharry Sickle Cell Specimen Collection.....	IV - 19
XIV.	Additional Information, Charts and Forms	
	Chart IV - 2, Pounds and Ounces to Grams Conversion.....	IV - 22
	Chart IV - 3, Military Time Conversion	IV - 23
	Chart IV - 4, Tennessee County Code List	IV - 24
	Newborn Screening Refusal Form.....	IV - 25
	Release of Records Request Form	IV - 26
	TDH Requisition for Laboratory Supplies Form	IV - 27

Newborn Screening
Laboratory (615) 262-6352
Results and Follow-Up (615) 262-6304

I. Introduction

Tennessee law requires that a blood sample shall be obtained from each infant born in the state, regardless of age, before discharge from the hospital and tested for specific genetic disorders. These genetic disorders can cause mental retardation or death if not treated quickly. To prevent the early affects of these disorders, the blood sample should be drawn in the infant's first 2 to 5 days of life. Infant's who are screened before they are 24 hours of age must be rescreened within 24 to 48 hours by a local health department, or private physician. Infant's who are born in a non-hospital setting must be taken to a hospital, local health department, or private physician between 24 and 48 hours of birth to have the blood sample collected. A few drops of blood from the infant's heel are absorbed into a special filter paper attached to the Newborn Screening Form PH-1582 and sent to the Newborn Screening Laboratory at the Tennessee Department of Health (TDH) Laboratory Services. The current protocol includes testing the blood sample for the following disorders:

PHENYLKETONURIA (PKU)

Infants with PKU lack an enzyme called phenylalanine hydrolase, which is needed to breakdown part of the protein in food. This protein part is called phenylalanine, and unless the condition is detected and treatment is initiated soon after birth, phenylalanine builds up in the blood and body tissues causing mental retardation. Infants are placed on a lifelong special diet low in phenylalanine, which can prevent mental retardation and other effects of PKU. Testing began in 1968 for this disorder, which has an incident rate of 1:12,000 infants.

CONGENITAL HYPOTHYROIDISM (TSH)

Hypothyroidism occurs when the body does not produce enough thyroid hormone from the thyroid gland. This hormone is called thyroxine, which is needed for brain and body growth. A decreased amount of thyroxine can interfere with normal growth and can lead to mental retardation. If detected early and hormone replacement is initiated, normal growth and development can take place. Testing began in 1980 for this disorder, which has an incident rate of 1:4,500 infants.

HEMOGLOBINOPATHIES

Hemoglobin is the part of red blood cells that carries oxygen. Hemoglobinopathies are diseases that affect the kind or amount of hemoglobin a person has in the red blood cells. Some hemoglobinopathies can cause anemias or thalassemias.

Infants of **Asian** or Mediterranean background, such as Italian, Greek, Arab, or Kurdish, should be rescreened for thalassemias at the 12-month routine child-health visit. These specimens can be sent to the Meharry Sickle Cell Center or other laboratory.

Newborn Screening (Continued)

Sickle cell disease is the most common hemoglobinopathy. The red blood cells are sticky and crescent or sickle-shaped and therefore do not move easily through the vascular system, decreasing the vital levels of oxygen carried throughout the body. There is a dramatic decrease in infant mortality and life-threatening complications when an infant is identified early and treated with antibiotics. Testing began in 1988 and the incidence rate is 1:350 in the Africa-American population. People of Hispanic, Asian, Arabic, or Mediterranean decent are also more likely to have hemoglobinopathies.

GALACTOSEMIA (GAL)

Infants with galactosemia lack an enzyme, uridyl transferase, needed to breakdown galactose, a kind of major sugar found in milk. Due to either an absent or low enzyme, galactose accumulates in the body leading to mental retardation, growth deficiency, blindness, overwhelming infection and death. Infants with galactosemia can rapidly become sick after only a few days of normal feeding. Infants who are detected early are placed on a lifelong special galactose-free diet. Testing began in 1992 for this disorder which has an incident rate of 1:60,000 infants.

CONGENITAL ADRENAL HYPERPLASIA (CAH)

In order for this test to be reported accurately it is imperative that the **CURRENT weight** of the infant be recorded correctly on the Newborn Screening Form. This disorder results from a deficiency in one or another of the enzymes of steroid syntheses. The deficiency results in an abnormal accumulation of hormones important for normal growth and development of puberty. Most affected infants have a second defect that affects normal electrolyte balance. Infants having this deficiency can become very sick within 1 to 2 weeks of birth. Female infants with this disorder may be assigned the wrong sex at birth. Male infants are often not identified until they are in crisis. With early diagnosis and treatment it is possible to achieve normal growth and development of puberty. Because this disorder is complex, infants may need to be monitored by a pediatric endocrinologist. Testing began in 2000 for this disorder and has an incident rate of 1:10,000 to 1:15,000 infants.

II. Test Methods

- Phenylketonuria and galactosemia tests are performed by enzymatic fluorometric methods. A manual fluorescent test is performed on all borderline and abnormal galactosemia results to detect the uridyl transferase enzyme activity. Enzyme activity is reported as normal, low, or absent.
- Thyroid testing is performed by a fluoroimmunoassay (FIA) method which detects the amount of thyroid stimulating hormone (TSH) present.
- Hemoglobin testing is performed by an isoelectric focusing (IEF) method.
- Congenital adrenal hyperplasia employs a fluoroimmunoassay (FIA) methodology, which detects the amount of 17 α -hydroxyprogesterone (17- α OHP).

III. Specimen Collection

A. Collection Form

Use the Newborn Screening Collection Form PH-1582. Forms are available from your local county health department. Health departments can order more forms from the Shipping Department of Laboratory Services. (See Page IV - 27 for the TDH REQUISITION FOR LABORATORY SUPPLIES FORM.) The expiration date of the filter paper is printed near the bottom right corner of the form in red ink. It will have the month/year. Forms are good until the last day of the month printed on the form for the year specified. (Example: EXP. 8/2000 means do not use the form after 8/31/00.) Blood collected on forms after the expiration date will be reported out as "*Unsatisfactory Filter Paper Expired*" and another specimen will have to be submitted.

B. Important Information

1. Infants more than 6 months old

The test methods used by the Newborn Screening Laboratory are *not suitable for infants greater than 6 months of age*. If the infant is greater than 6 months of age, contact the Metabolic Center closest to the provider to inquire as to what tests need to be performed and where to send the specimens.

2. Procedures When Infants Are Transfused Prior to the Newborn Screening *

If possible, collect a specimen for the newborn screen **before a transfusion** even if the infant is less than 24 hours of age. The hemoglobin will be accurate, and if it is normal, there is no need for follow-up hemoglobin testing.

The PKU, TSH, CAH, and Total GAL test methodologies and results (i.e. measurement of the metabolite) can be affected by transfusions. Any infant that was transfused within 72 hours prior to specimen collection with any of the following blood products will need to have their screen repeated:

- Exchange transfusion
- ECMO procedure
- Whole blood
- Plasma

GAL and Hemoglobinopathy results can be affected by red blood cell transfusions.

* If an infant has symptoms such as vomiting, diarrhea, dehydration, or jaundice, the newborn screen should be repeated immediately regardless of the number of days that have passed since the transfusion. In addition, the regional metabolic center or endocrinologist (depending on symptoms) should be contacted.

Repeat PKU, TSH, and CAH between 4 and 14 days after transfusion. **Repeat GAL** between 10 and 90 days after transfusion.

Collect a specimen for PKU, TSH, CAH, and GAL as soon as possible after the tenth day posttransfusion. Send the specimen to the TDH Newborn Screening Laboratory.

Newborn Screening (Continued)

Repeat Hemoglobinopathy at three months posttransfusion. Send the specimen to Meharry Sickle Cell Center. See Newborn Screening Section XIII for information about submitting specimens to Meharry Sickle Cell Center.

If the infant has been transfused, mark transfusion on the collection form and give the last transfusion date.

3. Infants in the NICU and Premies

Sick or premature newborns should have a specimen collected on or near the seventh day of age regardless of feeding status or before transfusion. PKU results are based upon the assumption the infant has had protein feed. Galactosemia results are based on the assumption the infant is on a lactose feeding. If the infant is sick or not feeding well and the physician feels a test is not accurate due to the infant's feed status, the physician is encouraged to obtain a repeat specimen. If an infant is on soy formula, the provider can request an enzyme to be done by writing "DO ENZ" next to the Gal test on the collection form.

Newborn Screening (Continued)

C. Instructions for Filling out the Newborn Screening Form PH-1582 (Revised 2/00)

It is important to fill in all of the blank lines on the form completely, legibly, and accurately. Use a ballpoint pen to legibly print the information on the form. If critical areas on the form are left blank the Maternal and Child Health Follow-up Section will have to call to obtain the information. This may cause a delay in reporting results.

FRONT

NEWBORN SCREENING		TO AVOID RECOLLECTION- Accurately complete the entire form	
Specimen	First Repeat <input type="checkbox"/>	Previous TDH#	Reason: () 24 hr. () Unsatisfactory () Abnormal () Transfused Test: () PKU () TSH () Hgb () Gal () CAH
Infant's Last Name	First	Previous Last Name	SEX: () 1. Male () 2. Female
Birth Date	Time of Birth	MIL TIME () 1. Single Birth () 2. Twin (A or B) () 3. Other	RACE: () 1. White () 4. Am. Ind. () 2. Black () 5. Other () 3. Asian () 6. Hispanic
Date Collected	Time Collected	Hospital of Birth	ETHNICITY: () 1. Hispanic () 2. Non-Hispanic
Phone	Spec. Collected By	Hospital Collected	MOTHER'S INFORMATION
Name	Address	Medical Record Number	Mother's Current Last Name
City	State	Lab	First
Zip Code	City	Lab	Age
	Zip Code	Lab	Address
		Lab	City
		Lab	State
		Lab	Zip Code
		Lab	Phone
		Lab	TennCare: () 1. Yes () 2. No
		Lab	Do 2 Dgt. Code Code
		Lab	City of State
		Lab	REASON
		Lab	Screening Screen
		Lab	W-101
		Lab	W-102
		Lab	W-103
		Lab	W-104
		Lab	W-105
		Lab	W-106
		Lab	W-107
		Lab	W-108
		Lab	W-109
		Lab	W-110
		Lab	W-111
		Lab	W-112
		Lab	W-113
		Lab	W-114
		Lab	W-115
		Lab	W-116
		Lab	W-117
		Lab	W-118
		Lab	W-119
		Lab	W-120
		Lab	W-121
		Lab	W-122
		Lab	W-123
		Lab	W-124
		Lab	W-125
		Lab	W-126
		Lab	W-127
		Lab	W-128
		Lab	W-129
		Lab	W-130
		Lab	W-131
		Lab	W-132
		Lab	W-133
		Lab	W-134
		Lab	W-135
		Lab	W-136
		Lab	W-137
		Lab	W-138
		Lab	W-139
		Lab	W-140
		Lab	W-141
		Lab	W-142
		Lab	W-143
		Lab	W-144
		Lab	W-145
		Lab	W-146
		Lab	W-147
		Lab	W-148
		Lab	W-149
		Lab	W-150
		Lab	W-151
		Lab	W-152
		Lab	W-153
		Lab	W-154
		Lab	W-155
		Lab	W-156
		Lab	W-157
		Lab	W-158
		Lab	W-159
		Lab	W-160
		Lab	W-161
		Lab	W-162
		Lab	W-163
		Lab	W-164
		Lab	W-165
		Lab	W-166
		Lab	W-167
		Lab	W-168
		Lab	W-169
		Lab	W-170
		Lab	W-171
		Lab	W-172
		Lab	W-173
		Lab	W-174
		Lab	W-175
		Lab	W-176
		Lab	W-177
		Lab	W-178
		Lab	W-179
		Lab	W-180
		Lab	W-181
		Lab	W-182
		Lab	W-183
		Lab	W-184
		Lab	W-185
		Lab	W-186
		Lab	W-187
		Lab	W-188
		Lab	W-189
		Lab	W-190
		Lab	W-191
		Lab	W-192
		Lab	W-193
		Lab	W-194
		Lab	W-195
		Lab	W-196
		Lab	W-197
		Lab	W-198
		Lab	W-199
		Lab	W-200

BACK

TO INSURE ACCURATE REPORTING AND FOLLOW-UP ALL BLANKS ON FORM MUST BE COMPLETED.
GALACTOSE RESULTS ARE BASED UPON THE ASSUMPTION THAT THE INFANT HAS HAD LACTOSE FEEDING.
UNLESS THE TRANSFUSION IS MARKED, THE ASSUMPTION IS THAT THE INFANT HAS NOT BEEN TRANSFUSED.

INSTRUCTIONS FOR SPECIMEN COLLECTIONS

1. Hold infant's limb in a dependent position to increase blood flow.
2. Clean heel thoroughly. Wipe with alcohol and dry before puncturing.
3. Puncture heel with sterile lancet deep enough to assure free flow of blood.
4. Wipe away first drop and discard.
5. Allow a large drop of blood to form on infant's heel. Apply the back side of the filter paper directly to the drop of blood at the puncture site (NOT the heel). **The drop of blood should be large enough to approximately fill one circle and must completely saturate through the paper. It should look the same on both sides of the filter card.**
DO NOT: a) Apply more than one drop of blood per circle.
DO NOT: b) Apply blood to both front and back of filter paper.
6. Apply blood to all circles.
7. Allow blood spots to completely dry in a horizontal position at room temperature (see diagram). Do not stack specimens while specimen is exposed. After drying, rewrap this cover sheet to its original position to protect specimen.
8. Mail within 24 hours of collection to:
NEWBORN SCREENING
DEPARTMENT OF HEALTH
LABORATORY SERVICES
P.O. BOX 309130
NASHVILLE, TENNESSEE 37230 - 5130

Newborn Screening (Continued)

1. Infant's Information

- a. Specimen First
Repeat

If this is the "first" specimen ever collected on the infant, place a mark in the "First" box. If the infant has previously had a newborn screen performed, at any time, then mark the "Repeat" box.

b. Previous TDH#

If this is a repeat test and the provider or the parent was sent a letter indicating that the specimen should be repeated, write the 11 digit unique TDH number included in the letter in this space. If possible enclose a copy of the letter with the repeat specimen.

c. Reason: () 24 hr. () Unsatisfactory () Abnormal () Transfused

If this is a repeat test and the provider or the parent was sent a letter indicating that the specimen should be repeated, a reason would be given on the letter.

If the first specimen was collected before the infant was 24 hours of age and this is why the repeat is being done, mark **(X) 24 hr.**

If the specimen was unsatisfactory for any reason and this is why the repeat is being done, mark **(X) Unsatisfactory.**

If the specimen had an abnormal test result and this is why the repeat is being done, mark **(X) Abnormal.**

If the specimen was collected after the infant was transfused and this is why the repeat is being done, mark **(X) Transfused.**

d. TEST: () PKU () TSH () HGB () GAL () CAH

If you marked the "Reason" as <24, place a mark next to the following test names: PKU, TSH, GAL, and CAH.

If the reason was Unsatisfactory, place a mark next to the following test names: PKU, TSH, HGB, GAL, and CAH.

If the reason was Abnormal, place a mark next to the name of the test that was reported as Abnormal.

If the reason was Transfused, place a mark next to the name of the test that was reported as transfused. See information regarding transfusions on Page IV - 3.

If the specimen is marked as a "First" specimen all test are automatically performed.

If the physician or healthcare provider has requested an Enzyme on the Galactose, write in the () GAL blank "DO ENZ". Example: (DO ENZ GAL).

Newborn Screening (Continued)

e. Infant's Last Name:

Legibly print the infant's last name.

f. Infant's First Name:

If the infant has a first name, print it legibly here. If the infant does not have a first name at the time of collection and it was a single birth, write "BOY" for the first name if it's a male infant or "GIRL" for the first name if it is a female infant. If there are multiple births also indicate the birth order by using A, B, C, etc. EXAMPLE: GIRL "A", BOY "B", GIRL "C", BOY "D".

g. Previous Last Name:

If the infant has had a change in their last name, legibly print their previous last name here. *This is very important!* Without this information we are unable to identify when an infant has had a repeat specimen collected.

h. Birth Date:

The day the infant was born. Write the date as MM / DD / YY. The date should be the same day as recorded on the infant's birth certificate.

i. Time of Birth:

The time the infant was born in military time. (See Page IV - 23 for a MILITARY TIME CONVERSION CHART.)
Example: Write 3:00 AM as 0300 and 3:00 PM as 1500. The use of strict military time will indicate AM or PM. The time should be the same time as recorded on the infant's birth certificate.

j. Date Collected:

This is the date the specimen was collected. Write the date as MM / DD / YY.

k. Time Collected:

The time the specimen was collected in military time. (See Page IV - 23 for a MILITARY TIME CONVERSION CHART.) Example: Write 6:00 AM as 0600 and 6:00 PM as 1800. The use of strict military time will indicate AM or PM.

l. Birth Status:

If the infant was a single birth, mark (X) 1. Single Birth. If the infants were twins mark (X) 2. Twin, and mark either () A for the first born or () B for the second born. If the delivery is triplets (or more), mark (X) 3. Other, followed by the number and letter to indicate the birth order. Example: Triplets would be written as 3A, 3B, and 3C. Quadruplets would be written as 4A, 4B, 4C, and 4D.

Newborn Screening (Continued)

m. Hospital of Birth and Hospital Collected:

THE REPORT WILL BE MAILED TO THE HOSPITAL OF COLLECTION.

Enter the seven-digit hospital code that indicates the location of the infant's birth and the collection hospital. The hospital code is available from the hospital of birth. This blank **must be completed** regardless of the provider. Private physicians and county health departments must also give the code for the hospital of birth on **each specimen submitted to the TDH Laboratory**. If the infant has been transferred, the hospital of collection may be different from the hospital of birth. **Make sure you have the correct codes for both hospitals.**

If the infant was not born in a Tennessee hospital, indicate the location by giving the name and location of the hospital of birth.

If the infant was not born in a hospital, enter the two-digit code for the county and record "HOME" in hospital block as shown. See Page IV-24 for a COUNTY CODE LIST FOR TENNESSEE. (For example Shelby County - Home Birth - 79 | HIOIME.)

n. Medical Record Number:

Give the unique patient number assigned to the infant in the hospital, doctor's office, or local health agency.

o. Sex:

If the infant is a boy mark (X) 1. Male. If the infant is a girl mark (X) 2. Female.

p. Race:

Place a mark next to the category which best reflects the race of the infant () 1. White, () 2. Black, () 3. Asian, () 4. American Indian, or () 5. Other.

q. Ethnicity:

Place a mark next to the category which best reflects the ethnicity of the infant () 1. Hispanic, or () 2. Nonhispanic.

r. Status of the infant at time of collection:

Transfused: If the infant was not transfused mark (X) No.

If the infant was transfused, in utero or after delivery, mark (X) Yes and then enter the **last date** of transfusion (MM / DD / YY) prior to the specimen being collected. The transfusion information must be recorded accurately so that the hemoglobinopathy, galactosemia, and other test results are accurate.

Newborn Screening (Continued)

Birth Weight: ***NOTE*** Give the weight of the infant at the time the specimen was collected. The current weight of the infant must be recorded accurately so the CAH test results are accurate. The weight must be recorded in grams. See Chart IV - 2 POUNDS AND OUNCES TO GRAMS CONVERSION.

Note: **Gestational age** is on older forms. If you are using one of these, mark through this and enter the current weight.

Feeding: Place a mark next to the method by which the infant is currently being fed () 1. Breast, () 2. Soy, () 3. Intra Venous, () 4. Lactose, () 5. Other.

2. Physician or Provider Information

The report will be mailed to the Physician or Provider listed here. The report will also be mailed to the hospital of collection.

a. Phone:

Legibly print the area code and phone number of the physician or health care provider to be contacted if there is an abnormal test result.

b. Spec. Collected By:

Legibly print the name of the individual collecting the specimen.

c. Name:

Legibly print the first and last name of the physician or health care provider.

d. Address, City, State, and Zip Code:

Legibly print the complete street address, city, state, and zip code + 4 of the physician or health care provider.

3. Mother's Information

Adoption: If a newborn has been adopted please write "ADOPTION CASE" on the form and put either the adoptive parents or the adoption agency's information in the spaces under Mother's Information. If a repeat specimen is required, letters will be sent to the mother listed on the form until the Newborn Screening lab receives a repeat specimen.

a. Mother's Current Last Name:

Legibly print the mother's legal last name.

b. First Name:

Legibly print the mother's legal first name.

Newborn Screening (Continued)

c. Age:

Legibly print the age of the mother at the time of the infant's birth.

d. Address, City, State, Zip Code:

Legibly print the mother's complete address where she is currently living including the city, state, and zip code + 4. If the mother's residence is not in the USA, name the country.

e. Phone Number:

Legibly print the area code and home phone number where the infant's mother can be contacted. If the mother does not have a phone, give the area code and phone number of a relative or neighbor who could easily contact the mother.

f. TennCare:

If the mother has TennCare, mark (X) Yes. If the mother has any other type of insurance mark (X) No.

g. Mother's Social Security Number:

Legibly print the mother's social security number.

h. County of Residence:

Legibly print the two-digit code that corresponds to the county in which the mother resides. (See Page IV - 24 for the COUNTY CODE LIST FOR TENNESSEE.)

4. Hearing Screen

Note: If you are only reporting hearing screens, **do not** use the newborn screening forms. Contact Jessica Trumble at 615-741-5234 for instructions regarding the correct method to report and for additional questions regarding the hearing screens.

Hospitals that provide newborn hearing screening utilizing the physiologic methods of ABR and/or OAE are requested to complete the Hearing Screen portion of the form. If your hospital does not provide newborn hearing screening, leave the spaces blank.

Place a mark next to ABR if the infant was tested by the Auditory Brainstem Response (ABR) or Automated Auditory Brainstem Response (AABR) method. Place a mark next to OAE if the infant was tested by the Otoacoustic Emissions (OAE) using the Distortion Product Otoacoustic Emissions (DPOAE) or the Transient Evoked Otoacoustic Emissions (TEOAE) method. If both methods were used, mark the method administered prior to discharge and the results.

Ear Screened: Place a mark next to R Ear if the right ear was screened. Place a mark next to L Ear if the left ear was screened. Mark both spaces if both ears were screened.

Newborn Screening (Continued)

Pass is the term used for a test that indicates the hearing is within normal limits. If the infant passed the test, place a mark next to _Pass. Refer is the term used for referral of infants for further evaluation. If the infant needs to be referred, place a mark next to _Refer.

5. Date Rec'd/Lab No (White area at bottom of form)

DO NOT WRITE IN THIS AREA

The laboratory will record the date and assign a specimen number to the specimen when it is received at the laboratory.


D. Filter Paper Collection Procedure

PRECAUTIONS:

Improper storage of the collection forms can result in the filter paper being altered by the absorption of ambient moisture in the atmosphere, or by compressing the filter paper if they are stacked on top of each other. Store the forms on their side in a cool, clean, dry place. Touching the filter paper, either before or after collection, with gloved or ungloved hands or substances such as powder, hand lotion, feeding formulas, or other materials can alter the filter paper and also the test results. Use of highlighters on the filter paper or wearing strong perfumes may also interfere with testing procedures.

Verify the patient information on the specimen collection form with the infant's identification band, and appropriate chart records before sticking the infant. Check the expiration date of the form. If it has expired do not use it. Appropriate biohazard precautions should be used. Wear powder free gloves and change gloves between infants. Never reuse lancets; place them in a puncture proof sharps container after use. Place biohazardous material in appropriate containers.

1. Select the puncture site on the heel. The shaded areas on the diagram indicate the preferred and least hazardous puncture sites. Draw an imaginary line between the fourth and fifth toes that runs parallel to the lateral aspect of the heel or a line running from the great toe that runs parallel to the medial aspect of the heel.

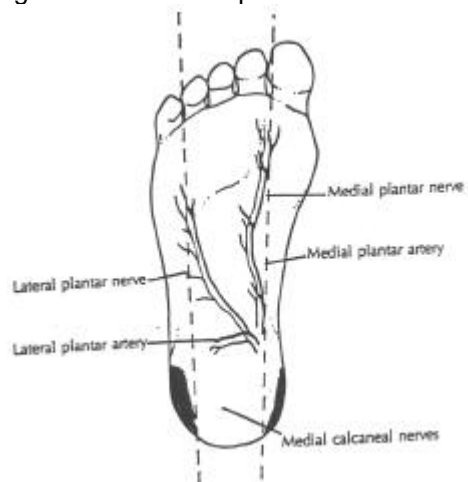

Newborn Screening (Continued)

2. Warming the heel with a warm (Temperature not to exceed 42 °C) wet towel or diaper or a commercially available heel warmer for not more than 3 minutes along with positioning the leg lower than the heart may dilate the blood vessels and therefore increase blood flow.
3. Disinfect the skin with 70% isopropyl alcohol and allow it to air dry. Vigorous rubbing may also stimulate blood flow to the area.
4. Puncture the skin with one continuous motion using a sterile automated lancet. On a full term infant the depth of the puncture should not exceed 2.0 mm. On premature or small infants the puncture should be less than 2.0 mm. Lancets with tips longer than 2.5 mm may cause excessive tissue damage, bone damage, or long term problems with walking.
5. Wipe away and discard the first drop of blood since it may be contaminated by alcohol or tissue fluid.
6. Allow the second drop (approximately 75-100 µl) to form by the spontaneous free flow of blood.
7. Touch the filter paper to the drop of blood as close to the center of the circle as possible. Allow the blood spot to expand within the circle and soak through to the other side. Turn the card over to make sure that the blood did soak through completely. DO NOT press or touch the filter paper to the puncture site.
 - a. Apply blood to only one side of the filter paper. It does not matter which side, however, the blood must soak through from one side to the other.
 - b. It is essential that only one drop of blood be used to fill a circle. If the circle cannot be filled with only one drop, go on to the next circle. Do not apply two drops to the same circle.
 - c. It is permissible for the blood to go outside the circles lines, but do not allow it to overlap on blood in an adjacent circle. Overfilling of the circle can cause supersaturation.
 - d. Do not use capillary tubes or syringes to fill circles. Such devices may alter the filter paper by scratching it, cause supersaturation, or blood clots.
 - e. Recollect the sample if tiny blood clots appear on the specimen or if any fluid or substance contaminates the specimen.
8. Once the blood collection is completed, hold the infant's foot above the heart level and press a sterile gauze to the puncture site until the bleeding has stopped. Adhesive bandages are not recommended.
9. Air dry the blood specimen horizontally at room temperature away from heat or direct sunlight for at least 3 hours. Do not allow the blood spots to touch any surface. Do not close the fold-over-paper protective flap for at least 3 hours and the specimen is completely dry.

IV. Shipment of Specimens

After the specimen has had time to dry for at least 3 hours, close the protective paper flap over the top of the blood spots. **DO NOT tape the flap closed or fold the form.** Place the form in a **paper** envelope. If you are sending more than one form, rotate the forms 180° so that the blood spots are not stacked directly over one another but are alternated. **DO NOT place the forms in any form of a plastic sealed bag.** This includes poly bags, ziplock bags, plasticene envelopes, or plastic shipping bags.

Mail or transport the specimen to the Tennessee Department of Health's Laboratory in Nashville within 24 hours of collection.

**If the specimens are going to be mailed
United States Postal Service Address:**

Newborn Screening
Tennessee Department of Health
Laboratory Services
P. O. Box 305130
Nashville, TN 37230-5130

**If the specimens are going to be
delivered by Courier, FedEx, or UPS:**

Newborn Screening
Tennessee Department of Health
Laboratory Services
630 Hart Lane
Nashville, TN 37216-2006

The U. S. Postal Service address is what appears on the back of the yellow fold-over flap. If you have your specimens sent by FedEx, UPS or another Delivery Company you must use the address on the right above.

V. Reporting Procedure and Interpretation

The results of all specimens are reported to the hospital that collected the specimen and the physician or provider listed on the form. If you need a copy of a Newborn Screening Report and your hospital or physician was not listed on the original form, you will need to submit a Release of Records request form signed by the legal parent or guardian. See Page IV - 26 for the REQUEST FOR RELEASE OF RECORDS FORM.

The laboratory reports all positive, or suspected positive, results to Maternal and Child Health, generally within 1 to 2 working days after the specimen is received in the laboratory. Maternal and Child Health then notifies the provider by telephone to initiate treatment of the patient, confirmation testing, and follow-up of the patient. The written report is mailed within 7 to 10 working days after receipt of the specimen.

Written reports of normal specimens are mailed within 7 to 10 working days after receipt in the laboratory.

All unsatisfactory specimens are tested, even though the integrity of the specimen is in question. If a positive, or suspected positive, is found, results are reported to Maternal and Child Health, so that the provider can be notified and treatment can be initiated. **However, all unsatisfactory specimens are reported as unsatisfactory and must be repeated.**

Newborn Screening (Continued)

Reference Ranges											
<p>Phenylketonuria (PKU) < 4.0 mg/dl Within Normal Limits 4.0 – 5.9 mg/dl Borderline 6.0 mg/dl or greater Positive</p>											
<p>Galactosemia (GAL) Positive Enzyme and a Total Gal < 15.0 mg/dl Within Normal Limits Positive Enzyme and a Total Gal ≥ 15.0 – 19.9 mg/dl Borderline Positive Enzyme and a Total Gal 20.0 mg/dl or greater Abnormal Low Enzyme and a Total Gal ≥ 7.0 - < 15.0 mg/dl BorLow Low Enzyme and a Total Gal ≥ 15.0 mg/dl AbLow Absent Enzyme and a Total Gal ≥ 5 mg/dl AbNeg</p> <p>It is assumed that the infant has had a lactose feeding and has not been transfused unless otherwise indicated on the Newborn Screening Collection Form.</p>											
<p>Congenital Hypothyroidism (TSH) <u>For infants 1 day through 7 days</u> < 33 μU/ml serum Within Normal Limits 33 – 55 μU/ml serum Borderline > 55 μU/ml serum Positive</p> <p><u>For infants 8 days through 6 months of age</u> < 13 μU/ml serum Within Normal Limits ≥ 13 μU/ml serum Positive</p> <p><u>For infants <24 hours of age</u> < 33 μU/ml serum Within Normal Limits but MUST REPEAT FILTER TEST ≥ 33 μU/ml serum Inconclusive Results MUST REPEAT FILTER TEST</p>											
<p>Hemoglobinopathy FA is a normal hemoglobin pattern for a young infant AF is a normal hemoglobin pattern for an old infant The result is reported as Within Normal Limits as long as no abnormal bands are detected.</p> <p>It is assumed that the infant has not been transfused unless otherwise indicated on the Newborn Screening Collection Form.</p>											
<p>Congenital Adrenal Hyperplasia (CAH)</p> <table style="margin-left: auto; margin-right: auto;"> <tr> <td>Birth Weight</td> <td></td> </tr> <tr> <td><1250 grams</td> <td><135 ng/mL</td> </tr> <tr> <td>1251-1750 grams</td> <td><90 ng/mL</td> </tr> <tr> <td>1751-2245 grams</td> <td><65 ng/mL</td> </tr> <tr> <td>>2250 grams</td> <td><50 ng/mL</td> </tr> </table>		Birth Weight		<1250 grams	<135 ng/mL	1251-1750 grams	<90 ng/mL	1751-2245 grams	<65 ng/mL	>2250 grams	<50 ng/mL
Birth Weight											
<1250 grams	<135 ng/mL										
1251-1750 grams	<90 ng/mL										
1751-2245 grams	<65 ng/mL										
>2250 grams	<50 ng/mL										

VII. Specimen Criteria

Newborn Screening Laboratory personnel closely examine each specimen that arrives in the lab for quality and quantity before testing begins. Two technologists look at all unsatisfactory specimens. The collection hospital, provider, and the infant's mother are notified by mail that a repeat specimen needs to be obtained. Unsatisfactory reports are usually mailed within 2 to 4 working days after the specimen is received. The specimens has to be recollected and sent to the laboratory which can be a costly delay if the infant is diagnosed with a genetic disorder. A satisfactory blood drop is one single drop of whole blood applied evenly and allowed to soak through the filter card and can be seen clearly with no white showing through on the other side. The spots must be large enough to allow FIVE 1/8-inch discs to be punched out with no white areas showing.

See Chart IV - 1 COMMON CAUSES OF UNSATISFACTORY NEWBORN SCREENING SPECIMENS for unsatisfactory specimen conditions and their possible causes.

Newborn Screening (Continued)

Chart IV - 1
Common Causes of Unsatisfactory Newborn Screening Specimens

Unsatisfactory Specimen Conditions	POSSIBLE CAUSE
Nonuniform	Applying many small drops of blood to each circle. Applying blood with any type of capillary tube. Touching the blood drops when they are wet. Uneven soaking through the filter card is caused by exposure to moisture, glove powder touching the filter card area before collection or the use of hand creams or lotions. Contaminated surfaces with any of the above contaminants may also cause nonuniform absorption of the blood.
Filter Paper Expired	Specimen was collected on filter paper that was expired.
> 6 months	The infant was greater than 6 months of age at the time the specimen was collected.
< 24 hours	The specimen was collected before the infant was 24 hours of age.
> 5 days	The specimen was received in the laboratory greater than 5 days after the date of collection. There may be a delay in the U.S. Postal service or a delay in the hospital mail service. Mail specimens within 24 hours of collection.
Contaminated	Specimen was contaminated in some way with something like alcohol, water, formula, urine, or hand lotion for example.
Inaccurate Information	Information on the form was inaccurate or incorrect. This most often occurs when the date or time of collection is written on the form as occurring before the date and time of birth.
Incomplete Information	All blanks on the form were not filled out completely.
Quantity Not Sufficient (QNS)	The drops of blood are too small. This can be caused by improper use of the lancet or dropping blood from a capillary tube device.
Supersaturated	The drops of blood are too large. The drops of blood overlap or touch one another. The filter card is pressed against the puncture site. The blood is dropped in very large drops from a capillary tube.
Cells & Serum Separated	The usual cause is squeezing the heel during the specimen collection. It can also be caused by waiting too long for the drop of blood to form or by clotted blood. When applying blood with a capillary tube device, if the blood is not well mixed it may appear on the filter card as clotted blood.
Clotted Specimen	Clotted specimens are due to improper puncture, application with a capillary tube device, and waiting too long for a drop of blood to form allowing the blood to clot.
Altered Card	Use of a capillary tube or syringe to apply the blood can scratch the filter card when wet, or rubbing the spot when it is still wet. It can also be caused by pressing the heel to the filter card during the collection process.
Both Sides	This is from applying blood to both sides of the filter card. It is easily recognized in the lab by holding the specimen up to a light and looking at both sides for shadows.
No Blood	The form was received but there was no blood collected on the filter paper.
Heated	Specimen appears much darker than usual and appears to have been heated. Caused by too long in transit especially during the summer. Heat and humidity can affect test results. Also caused by heating a specimen to dry it.
QNSCOM	Insufficient blood to complete testing. Quantity not sufficient for test completion.
Poly Bag	Specimen received in a poly bag, plastic zip lock bag, plasticene envelope, or plastic shipping bag

VIII. Confirmation Procedures

The laboratory reports a presumptive positive result to Maternal and Child Health (MCH) as soon as it is determined, generally within 24 to 48 hours after the specimen is received. MCH notifies the provider listed on the form by telephone and fax to initiate confirmatory testing, follow-up, and treatment of the infant. MCH also notifies the appropriate endocrinologist and Genetic Center or Sickle Cell Center. As soon as abnormal results are determined, MCH faxes or mails, by certified letter, presumptive positive results to the provider. Final results are mailed when all other test are completed, generally within 7 to 10 days after the specimen is received.

PKU Confirmation: Those infants who have an abnormal phenylalanine level on newborn screening require confirmatory testing for phenylketonuria (PKU). These infants should be referred immediately to the nearest genetic center as recommended by the Maternal and Child Health (MCH) Follow-up Program. The genetic center will provide information regarding follow-up procedures for blood collection and mailing.

Infants who have a borderline phenylalanine level are required to have follow-up. A repeat filter-paper blood specimen should be submitted to the TDH Laboratory.

Galactosemia Confirmation: Infants who require confirmatory diagnostic testing for galactosemia should be referred immediately to the nearest genetic center for instructions regarding blood specimen collection and mailing. The MCH Follow-up Program will provide information about the nearest genetic center.

Infants with a borderline result for galactosemia are required to have a repeat filter-paper blood specimen sent as soon as possible to the TDH Laboratory.

Infants who received a transfusion before the newborn screening test need to have a repeat filter-paper blood specimen collected ten days following the last date of transfusion and within three days of lactose feedings. Repeat galactosemia screening is done again three months after the transfusion.

Hemoglobinopathy Confirmation: All infants either transfused or whose initial screen identified an abnormal hemoglobin trait or disease must have a hemoglobin confirmation performed by the Meharry Sickle Cell Center. All supplies, including mailers, collection devices, and forms can be ordered from:

Comprehensive Sickle Cell Center
Meharry Medical College
Nashville, TN 37208

Telephone: (615) 327-6763

Congenital Hypothyroidism Confirmation: MCH will notify the provider of infants who require confirmatory diagnostic thyroid testing. MCH recommends serum thyroid testing for confirmatory testing. This may be handled by the primary care physician or an endocrinologist.

Congenital Adrenal Hyperplasia (CAH) Confirmation: Infants who require confirmatory diagnostic testing for CAH should be referred immediately to the nearest pediatric endocrinologist for instructions regarding blood specimen collection and follow-up.

Newborn Screening (Continued)

IX. Material for Parents

The state provides pamphlets for parents to educate them about the Newborn Screening Test. These yellow pamphlets can be obtained from Maternal and Child Health (MCH) by calling (615) 262-6304. The pamphlet answers questions most often asked and explains briefly the 5 disorders and how the testing is done.

There is also a green pamphlet "Your Baby Needs to be Rescreened" which is require by state law to be given to every parent whose infant leaves the hospital before 24 hours of age. You are still required to do a newborn screen on this infant before he/she leaves. On newborn screenings done before the infant is 24 hours old, only the Hemoglobin test is accurate, the other test will have to be repeated within 2 weeks. The green pamphlet explains this to the parents and also gives an overview of the 4 disorders the infant will need to be retested for.

X. Infant Discharged before Specimen Was Obtained

When the nursery or lab forgets to get a specimen on an infant before he/she leaves, we ask that you still send in a newborn screening collection form completely filled in and write, "LEFT BEFORE OBTAINED" on the form. You are responsible for notifying the parent and the physician of the need for the test. Newborn Screening will also contact the physician and make sure the infant is screened. When a collection form for these infants are not sent in, we are unaware that the infant exists.

XI. Refusal to Have Newborn Screening Performed

When a parent refuses the newborn screening test, we ask that you do send in a newborn screening collection form completely filled in and write "REFUSED" on the form, and we will follow up. If a parent refuses the test due to religious reasons, please fill out the NEWBORN SCREENING REFUSAL FORM (See Page IV - 25), have it notarized, and fax it to Maternal and Child Health at (615) 262-6458.

XII. Death Notice

If you are aware of an infant that has expired, please fax us the child's name, date of birth, mothers name, and the date the infant expired. If a specimen needs to be repeated we will continue to send letters to the mother until we receive a specimen. In this case, we do NOT want to notify the parent. Maternal and Child Health's fax number is (615) 262-6458.

Newborn Screening (Continued)

XIII. Specimen Collection for the Hemoglobin Electrophoresis Procedure from Meharry Sickle Cell Center

Proper collection and handling of specimens are essential when sending samples to the Meharry Sickle Cell Center. The following steps must be adhered to so specimens are not damaged in transit and comply with Occupational Safety and Health Administration (OSHA) standards and postal regulations. If these directions are not followed, your samples may not reach us. Your cooperation in this matter is appreciated.

Supplies Needed: Request form (Adult/Children), Microvette Tube EDTA, Plastic Biohazard Bags, Styrofoam Mailing Boxes, and Packing Tape.

Each Microvette tube contains enough EDTA to anticoagulant up to 200 µl of blood.

Note special instructions in tube package.


1. Label the Microvette tube with the patient's name (label not included).
2. Perform a heel stick using good laboratory practice and following the instructions of the lancet device manufacturer.
3. Wipe off the first drop of blood. Then touch the capillary tip of the Microvette Tube to the droplet of blood as it appears at the incision. Excessive pressure or squeezing on puncture area must be avoided. This could lead to premature coagulation in the capillary.
4. After the blood has been taken, the capillary unit is removed and properly discarded and the sample container is closed with the attached stopper.
5. Properly mix the blood with the EDTA. Hold the stoppered collection Microvette Tube between thumb and index finger and invert several times.
6. Complete the laboratory request form with all information. Include the address to send results. Fold the laboratory request form and place it in the unsealed pocket of a biohazard bag.
7. Place the sealed Microvette Tube in the ziplock pocket of the biohazard bag and place into a Styrofoam mailing box.
8. Place the Styrofoam box in the mailing sleeve and seal both ends. Please put the return address on the mailing box.
9. Mail appropriately packaged sample(s) to :

ATTN: Laboratory Supervisor
Comprehensive Sickle Cell Center
Meharry Medical College
1005 D. B. Todd Boulevard
Nashville, Tennessee 37208

Phone (615) 327-6763 Fax (615) 327-6008

Newborn Screening (Continued)

Microvette® CB 1000

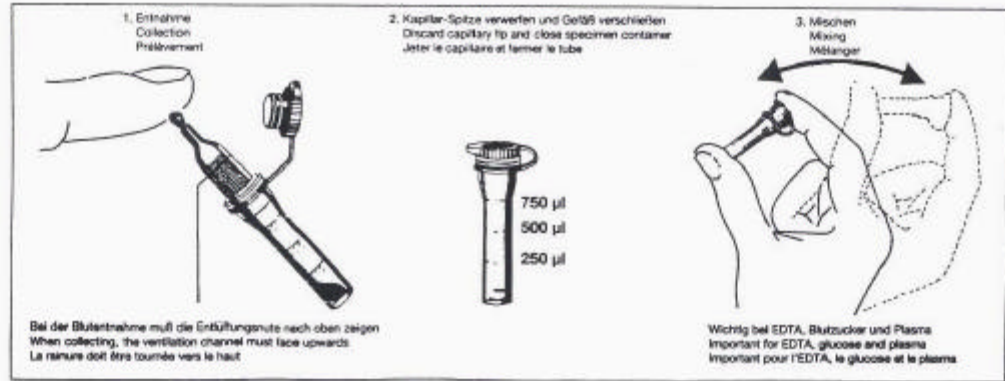

**Newborn Screening
Additional Information, Charts and Forms**

Newborn Screening (Continued)
 Chart IV - 2
 Pounds and Ounces to Grams Conversion

1 Kilo = 1000 grams

Example: To obtain grams equivalent to 5 pounds, 8 ounces, read "5" on the top scale, and read "8" on the side scale; it is equivalent to 2495 grams.

		POUNDS															
OUNCES		0	1	2	3	4	5	6	7	8	9	10	11	12	13	14	
	0	0	454	907	1361	1814	2268	2722	3175	3629	4082	4536	4990	5443	5897	6350	
	1	25	482	936	1389	1843	2295	2750	3203	3657	4111	4584	5018	5471	5925	6379	
	2	57	610	964	1417	1871	2325	2778	3232	3686	4139	4593	5046	5500	5953	6407	
	3	85	639	992	1448	1899	2353	2807	3260	3714	4157	4621	5075	5528	5982	6435	
	4	113	687	1021	1474	1928	2381	2835	3289	3742	4196	4649	5103	5557	6010	6464	
	5	142	695	1049	1503	1958	2410	2863	3317	3770	4224	4675	5131	5585	6038	6492	
	6	170	625	1077	1531	1964	2438	2882	3345	3799	4252	4706	5168	5613	6067	6520	
	7	198	652	1108	1559	2013	2466	2920	3374	3827	4281	4734	5188	5642	6095	6549	
	8	227	680	1134	1588	2041	2495	2948	3402	3856	4308	4763	5210	5670	6123	6577	
	9	255	709	1162	1616	2070	2523	2977	3430	3884	4337	4791	5245	5698	6152	6605	
	10	283	737	1191	1644	2098	2551	3005	3459	3912	4366	4819	5273	5727	6180	6631	
	11	312	765	1210	1673	2126	2580	3033	3487	3941	4394	4848	5301	5765	6209	6662	
	12	340	784	1247	1701	2155	2608	3062	3515	3969	4423	4876	5330	5783	6237	6690	
13	369	822	1278	1729	2183	2637	3090	3544	3997	4451	4904	5358	5812	6265	6719		
14	397	850	1304	1758	2211	2665	3118	3572	4026	4479	4933	5366	5840	6294	6747		
15	425	879	1332	1786	2240	2693	3147	3600	4054	4506	4961	5415	5868	6322	6776		

NOTE: 1 pound = 453.59237 grams;
 1 ounce = 28.349523 grams;
 1000 grams = 1 kilogram.

Gram equivalents have been rounded to whole number by adding one when the first decimal place is 5 or greater.

Chart IV - 3
Military Time Conversion Chart

CIVILIAN TIME		MILITARY TIME
1:00 AM		0100
2:00 AM		0200
3:00 AM		0300
4:00 AM		0400
5:00 AM		0500
6:00 AM		0600
7:00 AM		0700
8:00 AM		0800
9:00 AM		0900
10:00 AM		1000
11:00 AM		1100
12:00 a.m. noon		1200
1:00 PM		1300
2:00 PM		1400
3:00 PM		1500
4:00 PM		1600
5:00 PM		1700
6:00 PM		1800
7:00 PM		1900
8:00 PM		2000
9:00 PM		2100
10:00 PM		2200
11:00 PM		2300
12:00 p.m. midnight		2400

If the infant was born at 12:05 AM or 5 minutes after midnight, the military time would be written as 0005.

Newborn Screening (Continued)

Chart IV - 4
TENNESSEE COUNTY CODE LIST

COUNTY CODE	COUNTY NAME	COUNTY CODE	COUNTY NAME
01	ANDERSON	49	LAUDERDALE
02	BEDFORD	50	LAWRENCE
03	BENTON	51	LEWIS
04	BLEDSON	52	LINCOLN
05	BLOUNT	53	LOUDON
06	BRADLEY	54	MCMINN
07	CAMPBELL	55	MCNAIRY
08	CANNON	56	MACON
09	CARROLL	57	MASIDON
10	CARTER	58	MARION
11	CHEATHAM	59	MARSHALL
12	CHESTER	60	MAURY
13	CLAIBORNE	61	MEIGS
14	CLAY	62	MONROE
15	COCKE	63	MONTGOMERY
16	COFFEE	64	MOORE
17	CROCKETT	65	MORGAN
18	CUMBERLAND	66	OBION
19	DAVIDSON	67	OVERTON
20	DECATUR	68	PERRY
21	DEKALB	69	PICKETT
22	DICKSON	70	POLK
23	DYER	71	PUTNAM
24	FAYETTE	72	RHEA
25	FENTRESS	73	ROAN
26	FRANKLIN	74	ROBERTSON
27	GIBSON	75	RUTHERFORD
28	GILES	76	SCOTT
29	GRAINGER	77	SEQUATCHIE
30	GREENE	78	SEVIER
31	GRUNDY	79	SHELBY
32	HAMBLEN	80	SMITH
33	HAMILTON	81	STEWART
34	HANCOCK	82	SULLIVAN
35	HARDEMAN	83	SUMNER
36	HARDIN	84	TIPTON
37	HAWKINS	85	TROUSDALE
38	HAYWOOD	86	UNICOI
39	HENDERSON	87	UNION
40	HENRY	88	VAN BUREN
41	HICKMAN	89	WARREN
42	HOUSTON	90	WASHINGTON
43	HUMPHREYS	91	WAYNE
44	JACKSON	92	WEAKLY
45	JEFFERSON	93	WHITE
46	JOHNSON	94	WILLIAMSON
47	KNOX	95	WILSON
48	LAKE		

Newborn Screening (Continued)

Patient's Name _____
Mother's Name _____
Date of Birth _____
Hospital of Birth _____

NEWBORN SCREENING REFUSAL

I, _____ have been informed of the need for newborn screening which includes testing for phenylketonuria (PKU), congenital hypothyroidism, congenital adrenal hyperplasia, galactosemia, hemoglobinopathies, and other metabolic/genetic defects that would result in mental retardation or physical dysfunction as determined by the health department.

I have been informed that state law requires this testing and that violation of the law is a misdemeanor.

Nonetheless, I refuse these tests at this time for my newborn baby _____ because such test conflicts with my religious tenets and practices.

Under penalty of perjury, I affirm such refusal because of a conflict with my religious tenets and practices.

IN WITNESS WHEREOF, the undersigned has signed this form this _____ day of _____, 2_____.

Signature

STATE OF TENNESSEE

COUNTY OF _____

On this _____ day of _____, 2_____, before me personally appeared _____, to me known (or proved to me on the basis of satisfactory evidence) to be the person described in and who executed the foregoing instrument, and acknowledged that she executed the same as her own free act and deed.

Notary Public

My Commission Expires:

REFUSAL

Newborn Screening (Continued)

REQUEST FOR RELEASE OF RECORDS

Baby's Last Name _____

Baby's First Name _____

Date of Birth _____

Mother's Last Name _____

Mother's First Name _____

Mother's Social Security Number _____

Hospital of Birth or Collection _____

Any other previous last names used by the mother or the baby:

Name of Contact person in provider's office: _____

Direct Phone number of Contact Person: _(_____)_____

In the event we have trouble locating the baby listed above we will attempt to contact the provider for additional information.

I, _____ the legal parent or guardian of the baby listed above do hereby give consent for the State of Tennessee Newborn Screening Laboratory to release results to the Physician or Group Listed here _____.

Parent's Signature

Witness's Signature

**TENNESSEE DEPARTMENT OF HEALTH - LABORATORY SERVICES
REQUISITION FOR LABORATORY SUPPLIES**

SEND TO _____

ADDRESS _____

SIGNATURE _____ PHONE _____

*** **SUBMIT THIS FORM TO THE LABORATORY. BE SPECIFIC IN ORDERING.** ***

Quantity

COMPLETE KITS (Includes Form, Mailing Label, and Mailer)

REQUEST FORMS

	TB smear and culture	Miscellaneous Examination	PH-1573	
	Intestinal Parasites (O&P)	Water Bacteriology Analysis	PH-1575	
	Pinworm	Mycobacteriology Smear and Culture	PH-1577	
	Pertussis			
	Water Bacteriology	Syphilis Serology	PH-1578	
	Foodborne Outbreak Stool Collection Kit	Virology	PH-1579	
INDIVIDUAL ITEMS		Req. Laboratory Supplies	PH-1580A	
(Includes All Supplies Except Forms and Mailers)		Newborn Screening	PH-1582	
	Cary-Blair for Enteric/Bacterial Transport	(Order From Nashville Only)		
	Chlamydia/Gonorrhea Female Collection Kit	Gonococcus Culture	PH-1583	
	Chlamydia/Gonorrhea Male Collection Kit	Rabies	PH-1584	
	Gonorrhea Culture	Group A Streptococcus	PH-1587	
	Throat Cultures for Group A Strep (TSA)	Immunoserology	PH-1589	
	Vacutainer Tubes (7ml draw)	Rubella Serology	PH-1917	
	Viral Transport Media (1.5 ml. / tube)	HIV-1 Serology	PH-3173	
	Diphtheria (Pai) (Order From Nashville Only)	Chlamydia and Gonorrhea Detection	PH-3179	
	Charcoal Agar (For Pertussis Only)			

MAILING CONTAINERS

MAILING LABELS

	Medium (1-3 tubes, single bloods)	GONO-PAK	
	Large (3-6 tubes, multi-bloods, GC Culture)	MYCO	
		Specimen For Microbiological Examination	
		PH-0838	

SPECIAL NEEDS

ORDERING INFORMATION

THREE EASY WAYS TO ORDER

- ◆ Please order only what you will use in one month
- ◆ Please submit order one week in advance of need
- ◆ Be sure to check expiration dates carefully
- ◆ Please return all incomplete or expired materials

1. Mail to: Tennessee Department of Health
Laboratory Services
630 Hart Lane
Nashville, Tennessee 37247-0801
2. FAX (615) 262-6393
3. Call (615) 262-6322

FOR LABORATORY USE ONLY



RECEIVED _____ FILLED _____ INITIALS _____ SHIPPED _____
Date Date Date Date

PH-1580A
1160 REV 8/00

RDA

